

DRUGS OF ABUSE

APPROXIMATE* AMOUNT OF TIME DETECTABLE IN THE URINE

DRUG GROUP	ABBREVIATION	DETECTION PERIOD
AMPHETAMINE	AMP	1 - 3 DAYS
BARBITURATE	BAR	2 - 10 DAYS
BUPRENORPHINE	BUP	2 - 3 DAYS
BENZODIAZEPINE	BZO	1 - 6 DAYS +
CANNABIS	THC	14 - 84 DAYS / SINGLE USE 1-3 DAYS
COCAINE	COC	2 - 7 DAYS
METHADONE	MTD	1 - 4 DAYS
METHADONE METABOLITE	EDDP	1 - 4 DAYS
METHAMPHETAMINE	MET	3 - 5 DAYS
OPIATE	OPI	2 - 4 DAYS
PHENCYCLIDINE	PCP	7 - 14 DAYS
PROPOXYPHENE	PPX	1 - 3 DAYS


For further help or information call Quantum Diagnostics on 01992 651 111

**TO BE USED AS AN APPROXIMATE GUIDE ONLY*


READING RESULTS

**THERE IS NO MEANING ATTRIBUTED TO LINE COLOUR INTENSITY OR WIDTH
ANY EVIDENCE OF A LINE SHOULD BE CONSIDERED A LINE**


FOR DRUGS OF ABUSE URINE TEST STRIPS

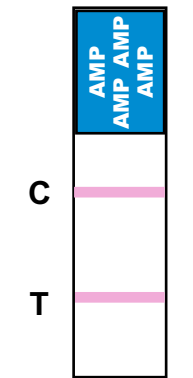


NEGATIVE

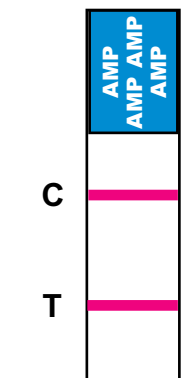


NEGATIVE

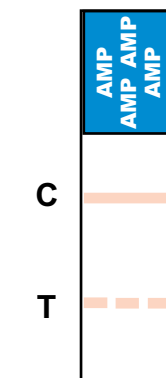


NEGATIVE

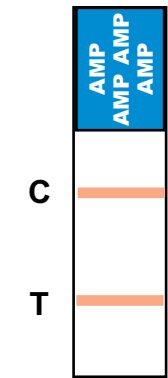


NEGATIVE

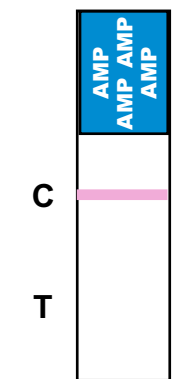


POSITIVE

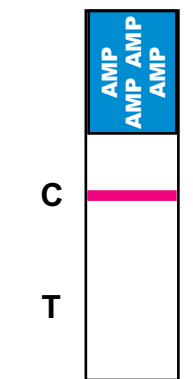

POSITIVE

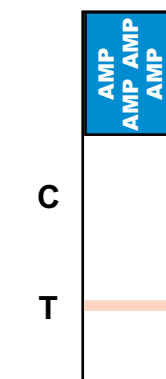

INVALID/RETEST

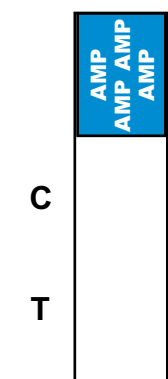

INVALID/RETEST